





OXYGEN INSTRUCTIONS

The on/off for the tank is on the top of the stem-----

A new tank will have a seal over the oxygen outlet-----


Under the outlet are two pin holes-----

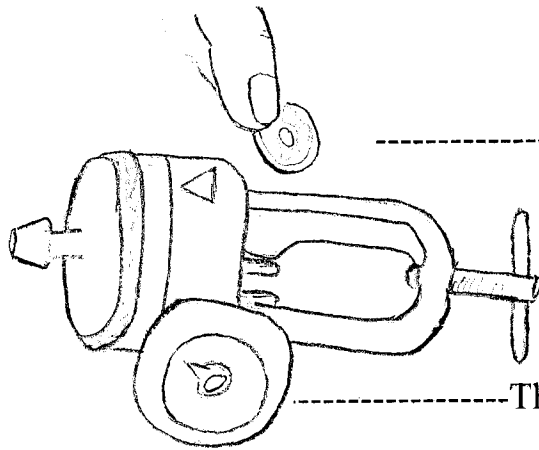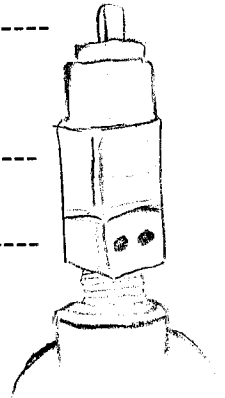


-----The regulator must have a washer
between it and the stem outlet.

-----The amount of oxygen in the tank is shown
on this pressure gauge.

Slide the regulator over the stem of the tank,-----
matching up the pins and the pinholes.
Snug the t-screw tightly.

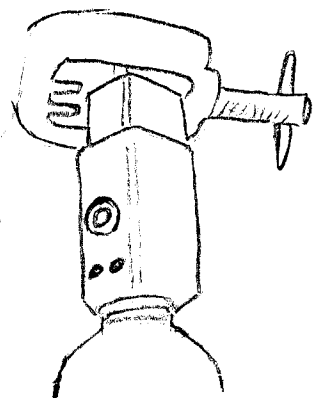


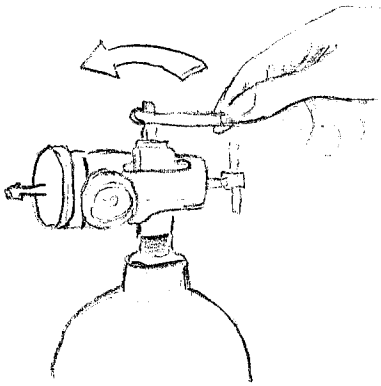

-----Using the wrench, turn the stem of the tank counterclockwise, to turn the tank on.

Set the flow to the prescribed flow rate-----

Connect the cannula to the outlet nipple-----

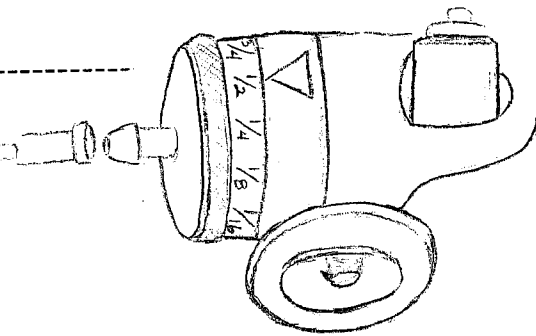

SAFETY PRECAUTIONS

Store cylinders lying down, and away from sources of heat and flame.

Avoid use of petroleum based products around oxygen.

Do not store oxygen cylinders in an unventilated area.

Do not smoke around your baby—or the oxygen equipment.

Check your oxygen supply often—please allow 24 hour notice for deliveries.

Always keep one backup cylinder for emergency use.

**CALL CHERUB MEDICAL SUPPLY
IF YOU HAVE PROBLEMS OR QUESTIONS**